
NS 5L SX

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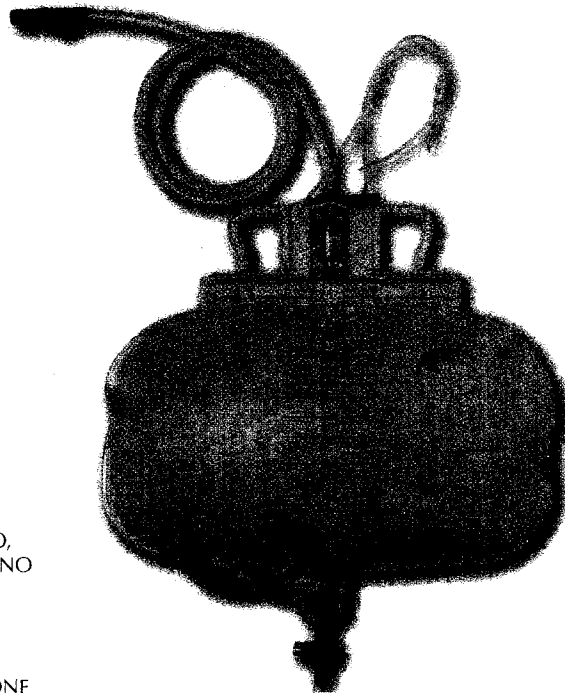
5000 cc
STRENGTHENED HANGER
PUNTO DI PRELIEVO
DRIP CHAMBER,
WITH NON RETOURNE VALVE
SLIDE CLAMP
VENT
T VALVE OUTLET
BLOCKING CLAMP
STRINGS FOR ATTACHING BAG

CODICE / REF.	DENOMINAZIONE / NAME	CARATTERISTICHE	DESCRIPTION	CONF. / PACKAGE
3808	NS 110 5L SX	TUBO 110 CM	TUBO 110 CM Long	50 Pz/Pcs
3809	NS 130 5L SX	TUBO 130 CM	TUBO 130 CM Long	

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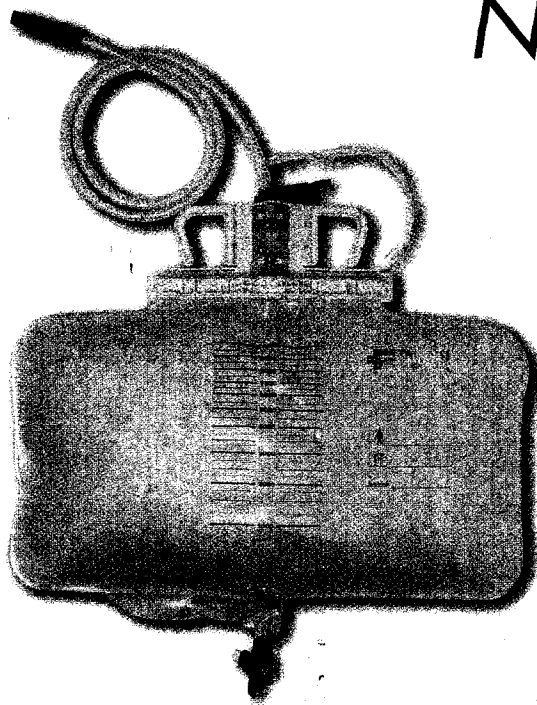10.000 cc
STRENGTHENED HANGER
PUNTO DI PRELIEVO
DRIP CHAMBER,
WITH NON RETOURNE VALVE
SLIDE CLAMP
SHEET CLIP
VENT
T VALVE OUTLET
BLOCKING CLAMP
STRINGS FOR ATTACHING BAG

CODICE / REF.	DENOMINAZIONE / NAME	CARATTERISTICHE	DESCRIPTION	CONF. / PACKAGE
3805	NS 110 10L SX	TUBO 110 CM	TUBO 110 CM Long	50 Pz/Pcs
3806	NS 130 10L SX	TUBO 110 CM	TUBO 110 CM Long	

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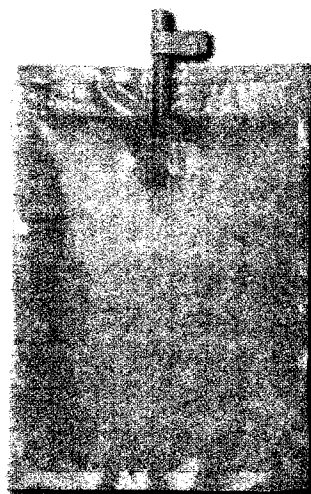
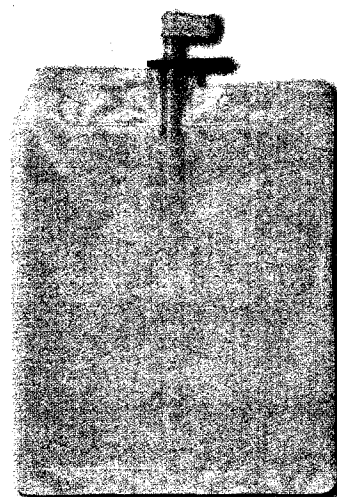
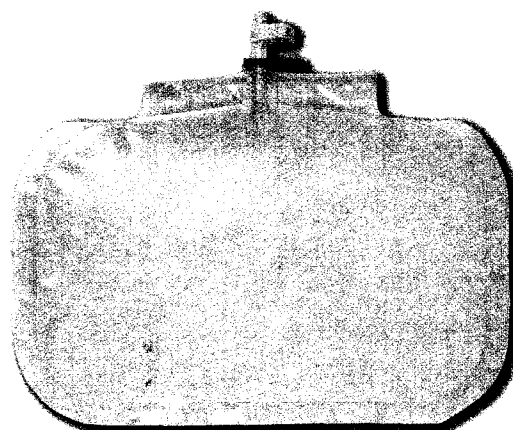
Meta
reserves the
right to modify
the products at
any time
without notice.

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
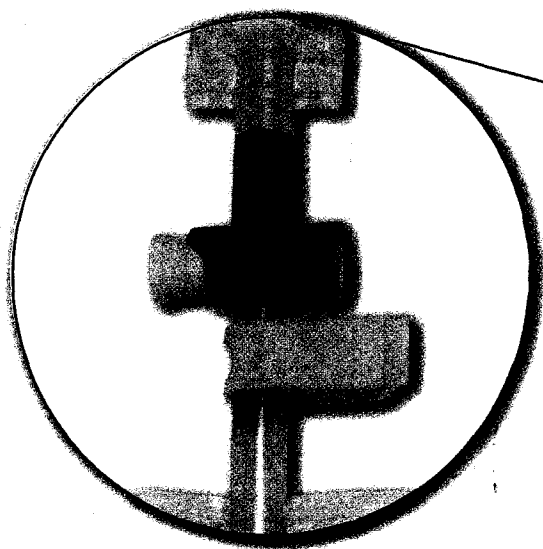
The New system Line (NS) identifies urine collecting bags and bags for vesical drainage endowed with a special system for emptying which contemplates the use of secondary emptying bags. In many cases, in fact, the emptying of bags still takes place inside open containers, which causes remarkable inconveniences and above all seriously compromises the hygiene of the ward favouring the transmission of infections due to the contact with urine. The New System Line easily solves these problems: by using the emptying bag Unibag to collect the urine coming from the closed circuit primary bag any contact with urine is avoided. The connection between the two bags is through a stopcock with a blocking clamp which guarantees a perfect seal.

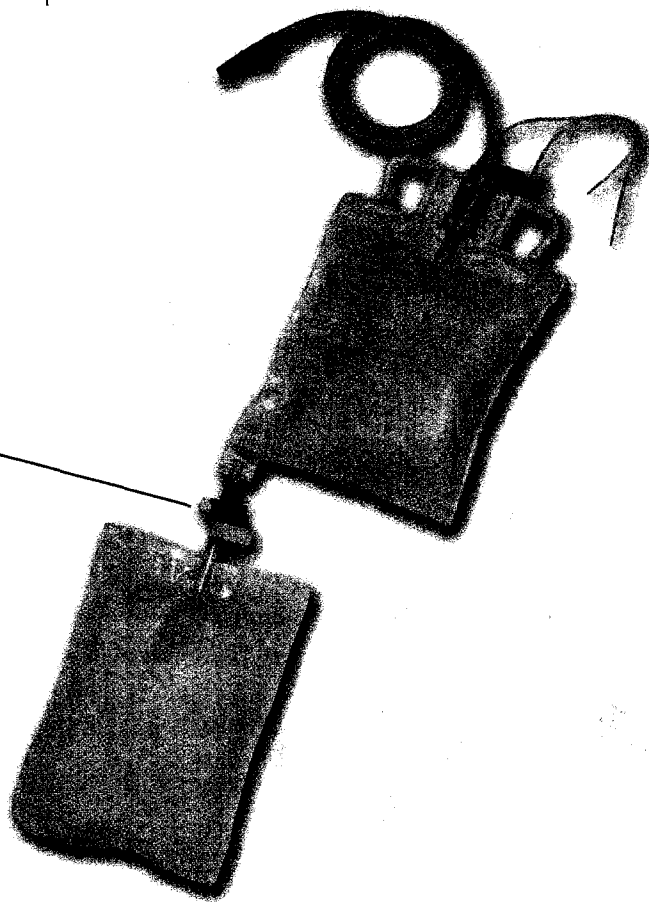
The New System Line are endowed with a special non-return valve. Usually urine bags have non-return valves whose functioning is based on a kind of gluing due to the contact with urine or on the displacement, due to the pressure of the urine going back, of rubber diaphragms.

In both cases the valve is to be considered as open to the going back up of bacteria since its being closed is determined by a likely counter pressure exercised by urine if the bag is bent or because of accidental knocks.

With these new series of bags also the problem of the compensation between drip chamber and bag is solved since the ventilation guaranteed through a unique filter allows the keeping of the same pressure in both areas.

Sacche Urina a Circuito Chiuso / Closed Circuit Urine Bags

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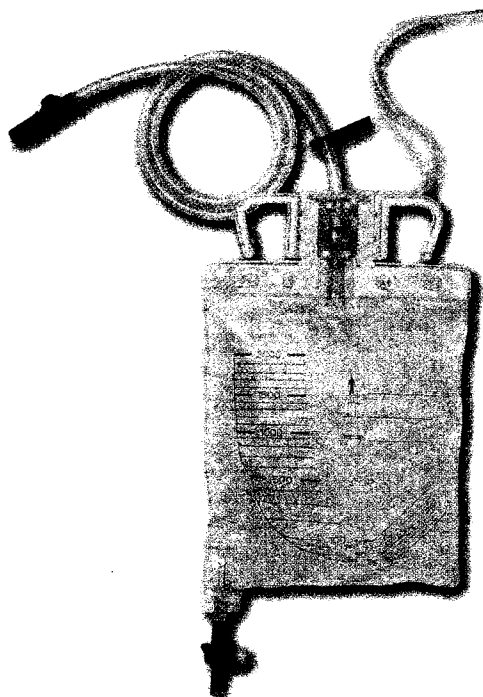
2000 cc
HANGER
SAMPLING SITE
NON RETOURNE VALVE
SLIDE CLAMP
VENT
T VALVE OUTLET
BLOCKING CLAMP
STRINGS FOR ATTACHING BAG

CODICE / REF.	DENOMINAZIONE / NAME	CARATTERISTICHE	DESCRIPTION	CONF. / PACKAGE
3813	NS 90 2L ST	TUBO 90 CM	TUBO 90 CM Long	55 Pz/Pcs
3814	NS 110 2L ST	TUBO 110 CM	TUBO 110 CM Long	
3815	NS 130 2L ST	TUBO 130 CM	TUBO 130 CM Long	

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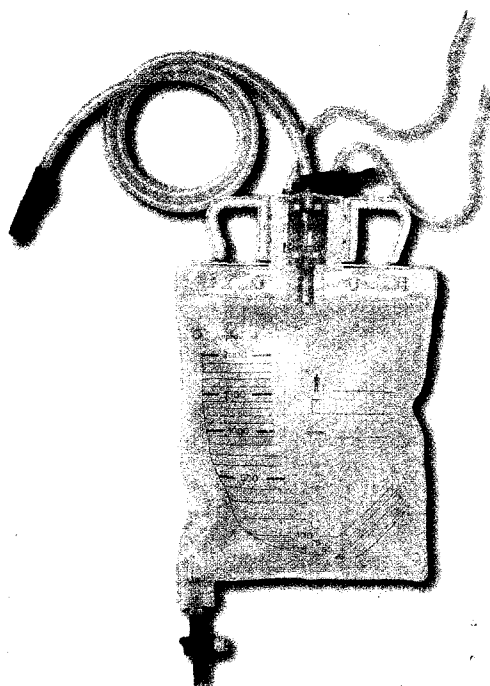2000 cc
HANGER
PUNTO DI PRELIEVO
DRIP CHAMBER,
WITH NON RETOURNE VALVE
SLIDE CLAMP
SHEET CLIP
VENT
T VALVE OUTLET
BLOCKING CLAMP
STRINGS FOR ATTACHING BAG

CODICE / REF.	DENOMINAZIONE / NAME	CARATTERISTICHE	DESCRIPTION	CONF. / PACKAGE
3810	NS 90 2L SX	TUBO 90 CM	TUBO 90 CM Long	55 Pz/Pcs
3811	NS 110 2L SX	TUBO 110 CM	TUBO 110 CM Long	
3812	NS 130 2L SX	TUBO 130 CM	TUBO 130 CM Long	

C.G.M. S.p.A.
Divisione
Medicale Meta
Via E. Villa, 7
42100
Reggio Emilia
Italy
Tel: +39 0522502311
Fax: +39 0522502333
www.metahosp.com
sales@metahosp.it